




Configure RASWin license for servers and clients:

Prerequisites: Sentinel ACC is required, this is installed either manually or automatically after RASWin installation.

Server side configuration


- 1) Start Sentinel Admin Control panel on the server using your web browser with this url: **http://localhost:1947.**
- 2) Go to **Configuration** on the left menu.
- 3) Click on tab **Users**.

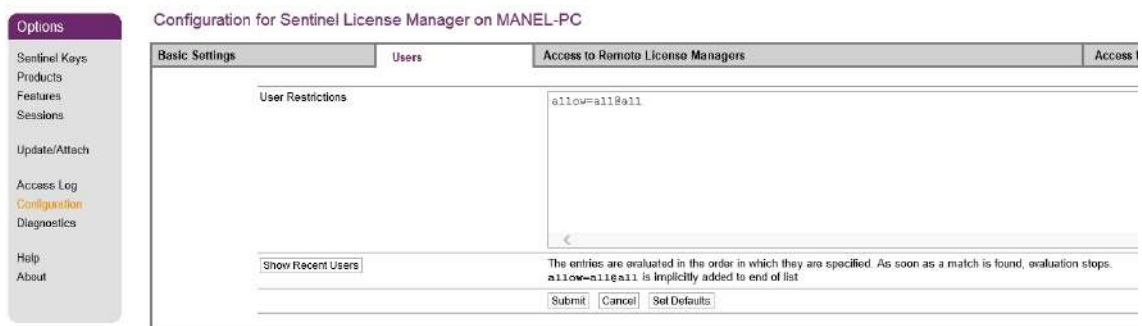

- 4) Ensure **User Restrictions** are allow=all@all
- 5) If you want to personalize access click on **Show Recent users**
- 6) On the list click each user you want to grant access by pressing “allow” or “deny”.

Client Side Configuration

- 1) On your Client machine start Sentinel Admin Control Panel opening the following URL: <http://localhost:1947>
- 2) Go to **Configuration** on the left menu.
- 3) Click on tab **Access To Remote License Managers**.
- 4) In the **Remote License Search parameters** enter the server ip.

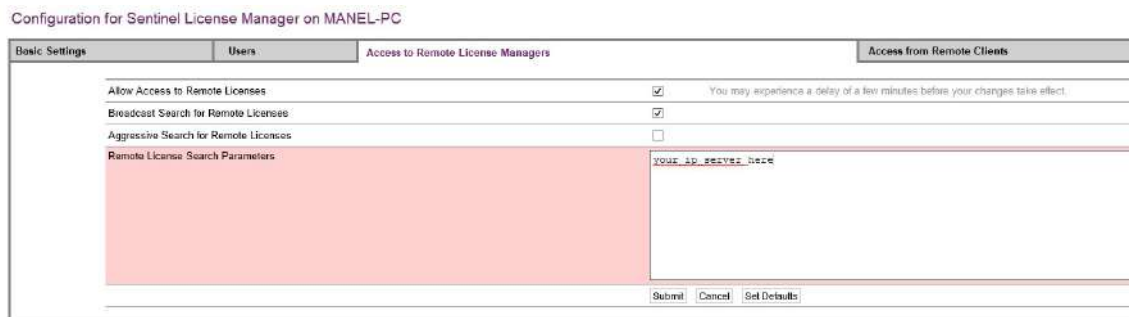

C:\Program Files (x86)\Common Files\AS\bin\Share\HAS\Phasplm.ini